
**CEN TEK®**

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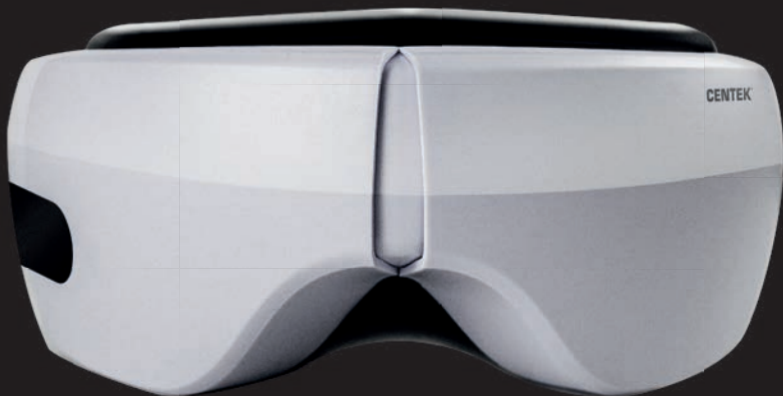# ENGLISH

## DEAR CUSTOMER!

Eye massage has contraindications. Consult a doctor before use.

The use of the eye massager is contraindicated:

- during pregnancy;
- if there is any skin damage in the area to be massaged, including burns and skin diseases;
- along with painkillers and sedative medications;
- for people using heart pacemakers.

Prior to operating the appliance, please carefully read the present manual which contains important information regarding its proper and safe operation. Take care not to lose the present manual. The manufacturer will not be responsible for any damage resulting from improper use of the appliance or from non-observance of the instructions and provisions of the present manual, as well as from attempts by unskilled persons at repairing the appliance. Should You decide to transfer the appliance to another person, please pass the present manual over along with it.

**The eye massager is a home appliance and is not designed for industrial applications!**

### 1. SAFETY PRECAUTIONS

The appliance purchased by You complies with all official safety standards applicable to electric appliances in the Russian Federation.

The following requirements should be strictly observed in order to avoid situations hazardous for human life and health, as well as premature failures of the appliance:

- Before connecting the appliance to the power supply grid, ensure that its specifications correspond to the values indicated on the appliance.
- Do not use the appliance if it is damaged or malfunctioning. In such cases, do not attempt to repair the appliance on your own. Apply to a service centre authorized by the Manufacturer.
- In case of power (USB) cable damage it should be replaced by the Manufacturer, service facility or similar skilled personnel.
- Do not use any accessories other than those recommended by the Manufacturer.
- Do not submerge the appliance in water or other liquids as this may damage the appliance.
- Handle the appliance with care, do not drop or expose it to strong mechanical impacts.
- Ensure the unit is turned off when not in use. Keep the appliance out of reach of children.
- Do not use the appliance while driving any kind of vehicle.
- Do not use the appliance while feeling unwell or sleepy.
- If You experience unpleasant sensations or indisposition during massage, stop the session and consult a doctor.
- Use the appliance only for its intended purpose; follow the instructions of the present manual.
- Do not use the massager outdoors and in indoor environments where any substances are sprayed out or rooms with high humidity.
- Do not cover the unit during operation as this may cause overheating and ignition. Protect the unit from direct sunlight, high temperatures and humidity.
- The appliance should not be used by physically, sensory or mentally handicapped persons (including children) or by persons without sufficient experience or knowledge, unless such persons are under supervision, or have been instructed on the use of the appliance, by the persons responsible for their safety. Take precautions not to children play with the appliance.
- The appliance has a heated surface. Persons insensitive to heat should be careful when

using the appliance.

- Special instructions for transportation (shipment) and sale: none.
- In case of a malfunction please apply to a service centre.

### CONTRAINDICATIONS:

- Severe chronic diseases in acute state
- Hypotonic conditions
- Infectious diseases
- Acute eye inflammatory diseases, glaucoma
- A doctor's advice is required in case of neurological diseases

### ATTENTION!

- Children under 8 years of age should only use the massager in the presence and under constant supervision of adults.
- Assume a convenient recumbent position prior to a massage session.
- Do not use the massager when walking or driving a car.
- If contact lenses are worn, remove them prior to using the massager.
- The duration of a massage session should not exceed 15 minutes.
- Do not use the massager while sleeping.
- Do not use the massager in bathrooms or other environments with high humidity.
- Use the massager only for its intended purpose.
- Do not disassemble the appliance yourselves.
- The device must only be powered by the voltage value marked on the device.

### 2. APPLIANCE DESCRIPTION

1. Massager
2. Power button
3. Fixing strap with buckle
4. Control panel
5. USB charger port
6. ON button
7. Massage mode indicator
8. Music Therapy mode indicator
9. Thermal Therapy mode indicator
10. Vibration Massage mode indicator
11. Battery charging indicator

### 3. SCOPE OF SUPPLY

- Massager - 1 pc.
- USB cable - 1 pc.
- Storage bag - 1 pc.
- Operating Manual - 1 pc.

### 4. OPERATING PROCEDURE

1. To turn the unit on press and hold the button for 3 seconds.
2. Select the desired modes by shortly clicking the button:

- Massage mode: massage + music
  - Thermal mode: massage + music + heating
  - Vibration mode: massage + music + heating + vibration
  - Sleep mode: music + heating
3. Bluetooth connection. Double-clicking the button activates the Bluetooth connection. Select **CENTER RELAX** from the device list and wait for the synchronization (no password is required).
  4. Для выключения прибора нажмите и удерживайте кнопку в течение 3 секунд.
  5. The massager is powered by a lithium battery. If you notice that the unit's output has significantly decreased, charge the massager.

#### RECOMMENDATIONS FOR USE

The use of the massager is recommended to:

Students who often suffer from pseudomyopia related to accommodation cramps caused by overfatigue.

Office workers whose work is related to high visual stress: PC operators, accountants, designers, secretaries, teachers, drivers, etc.

Persons suffering from presbyopia (age-related alteration of eye refracting power that require wearing positive diopter spectacles).

Persons with wrinkles, bags and dark circles around the eyes.

Persons suffering from headaches or eye pain, neuralgia or insomnia.

#### 5. MAINTENANCE

Wipe the massager body and control panel with soft flannel cloth. In case of bad contamination dampen the cloth in a soap solution. The use of solvents is not allowed! Protect the appliance from the ingress of water and other liquids. Store the appliance in a location protected from direct sunlight and dust. Do not store the appliance near heat sources and in places with high humidity.

Keep the appliance out of reach of children!

#### 6. SPECIFICATIONS

The appliance is powered by an internal battery. The battery is charged from a separate power supply unit with an output of 5 V  1.4 A connected to the appliance by a USB Type A interface.

- Power supply: 900 mAh Li-ion battery
- Output: 7 W
- Input voltage / current: 5 V  1.4 V
- Full charging time - 3 hours
- Battery life - up to 1 week on a full charge.

#### 7. DISPOSAL

Upon expiry of its service life the appliance can be disposed of separately from ordinary domestic waste. You can turn it over to a specialized electric and electronic appliance disposal and recycling centre.

#### 8. CERTIFICATION DATA, GUARANTEE

The service life of this item is 5 years since the date of sale to the end user, provided the item will be used in strict accordance with the present manual and the applicable technical

standards. Upon expiry of the service life please contact the nearest authorized service centre for suggestions as to further operation of the appliance. The manufacturing date of the appliance is integrated into the serial number (digits 2 and 3 - year, digits 4 and 5 - month of manufacture). Guarantee and post-guarantee maintenance is provided by the General Service Center OOO «Larina-Servis» (Russia, Krasnodar), phone: +7 (861) 991-05-42. The address of the facility authorized with receiving buyer claims, maintenance and repairs is shown on the website <https://centek.ru/servis>, otherwise refer to the store or other legal entity, where you have purchased this appliances.

The product is provided with the Certificate of Conformity.

#### 9. MANUFACTURER AND IMPORTER INFORMATION

Imported to EAEO by: Individualnyi Predprinimatel Asumryan Karo Shurikovich. Address: Russia, 350912 Krasnodar, ul. Atamana Lysenko, 23, pgt. Pashkovskiy, phone: +7 (861) 2-600-900.

#### DEAR VALUED CONSUMER!

The guarantee period for all appliances amounts to 12 months since the date of sale to the end consumer. With this Warranty Card the Manufacturer confirms the proper condition of the appliance and undertakes to repair, free of charge, all defects caused by the Manufacturer's default.

#### GUARANTEE MAINTENANCE TERMS

1. The guarantee shall be valid if the following documentation conditions are observed:

- The original Warranty Card shall be filled out correctly and clearly, indicating the of appliance model, serial number and sale date, the stamp of the seller company and signature of its representative shall be available, as well as the stamps on each coupon and the buyer's signature in the Warranty Card.

The Manufacturer reserves the right to reject guarantee maintenance if the buyer fails to provide the abovementioned documents or if the information contained therein is incomplete, illegible or ambiguous.

2. The guarantee shall be valid if the following operation conditions are observed:

- the appliance is used in strict accordance with the present manual;
  - the safety rules and requirements are observed.
3. The guarantee does not apply to regular maintenance, cleaning, installation and setup of the appliance at the owner's location.
  4. The guarantee does not apply in the following cases:
    - mechanical damage;
    - normal wear of the appliance;
    - non-observance of operation requirements or erroneous actions on the part of the user;
    - improper installation or transportation;
    - natural disasters (lightning, fire, flood, etc.), as well as any other causes beyond the Manufacturer's or Seller's control;
    - ingress of foreign objects, liquids or insects into the inside of the appliance;
    - repairs or modifications to the appliance by unauthorized persons;

- use of the appliance for commercial purposes (loads exceeding normal domestic application), connection of the appliance to utility and telecommunication networks incompliant with national technical standards;
  - failure of the following accessories, if replacement of these is allowed by design and does not involve disassembly of the appliance:
    - a) remote control units, rechargeable batteries, replaceable batteries (power cells), external power supply units and chargers;
    - b) consumables and accessories (packaging, covers, slings, carry bags, mesh screens, blades, flasks, plates, supports, grilles, spits, hoses, tubes, brushes, nozzles, dust collectors, filters, smell absorbers);
  - for battery-powered appliances - operation with inappropriate or depleted batteries;
  - for rechargeable battery-powered appliances - any damage caused by improper battery charging and recharging.
5. The present guarantee is provided by the Manufacturer in addition to the consumer rights determined by the applicable legislation and does not limit these rights in any way.
6. The Manufacturer shall not be made liable for possible damage, directly or indirectly inflicted by CENTEK brand products on people, domestic animals, the consumer's or third persons' property if such damage is inflicted as the result of non-observance of appliance installation and operation requirements, deliberate or negligent actions (omission) by the consumer or third persons, as well as of force majeure circumstances.
7. Service centers will only accept appliances in a clean state (there should be no food residues, dust and other contamination on the appliance).

**The Manufacturer reserves the right to change the design and specifications of the appliance without prior notice.**

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Attention! Please require the Seller to completely fill in the Warranty card and tear-off coupons. The Warranty Card is only valid if  
completed with all stamps and records.

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These tear-off coupons shall be filled in by the representative of the sales company

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